



SUNRISE PHARMA EQUIPMENTS

**MANUFACTURER & EXPORTER OF
PHARMACEUTICALS MACHINERIES**



**PROCESS & PACKAGING
FOR LIQUID & OINTMENT**

sales.sunrisepharma@gmail.com
www.sunrisepharma.in

SYRUP MANUFACTURING PLANT

We provide customized, end to end solutions for liquid oral plants that are fully automated with PLC controls or semi-automated as required. Our plants are designed with the highest quality of fabrication standards using the best components and instruments. They are also fully optimized for delivering the best performance with reduced life-cycle cost.



KEY FEATURES, OUR PRODUCTS AND SERVICES

- Design of process, equipment's and transfer piping
- Transfer systems for sugar through vacuum system
- Sugar preparation system with jacketed, insulated and cladded vessel with bottom entry propeller mixer for vigorous mixing
- In-line high shear mixer for suspensions
- Process vessels rated for vacuum
- Lobe pumps
- Electrically heat traced, insulated transfer lines of liquid glucose
- Back-washable stainless steel sintered filter for filtration
- Sanitary piping systems with butterfly valves
- Stainless steel operating platform with staircase and railing
- Compliant PLC based automation system
- Automated CIP system
- Highly experienced team for installation and validation support
- Capacity 500 -10,000 Ltrs

PROCESS EQUIPMENTS

- | | |
|----------------------------------|----------------------|
| • Material handling and transfer | • Platform structure |
| • Sugar syrup preparation | • Transfer panels |
| • Sugar syrup pre-filtration | • CIP Systems |
| • Solution preparation | • Plant automation |
| | • Solution storage |

OINTMENT MANUFACTURING PLANT



KEY FEATURES, OUR PRODUCTS AND SERVICES

- All contact parts are SS 316 L, and non-contact parts are SS 304.
- ointment Manufacturing vessels have an opening arrangement of top dish hydraulically operated.
- We use silicone food grade gaskets.
- We make the top-mounted mixer with SS 316 L anchor with a mirror polish, and PTFE arranged baffles/scrapper.
- The high-speed homogenizer at the bottom arrangement with a mechanical seal.
- High sanitary standard Design & Manufacturing.
- All material transfers are done by vacuum or by transfer pumps.
- The entire plant is equipped with CIP & SIP connections.
- The jacket has heating & cooling with a temperature controller.
- VFD for Anchor agitators to vary the RPM as per Recipe and required process.
- The vacuum system helps transfer product from one vessel to another vessel and zero retention in a pre-phase vessel & interconnecting pipeline.

PROCESS EQUIPMENTS

Wax vessel:

- Open top with top Entry agitator Cowl disc type or paddle with Slow Speed Anchor.
- Close top with top/bottom entry agitator Cowl disc type or propeller type stirrer.

Water Vessel:

- Open top with top entry propeller stirrer
- Close top with bottom or full entry propeller


Main Vessel:

- Choice of various types of anchor design.
- Homogenizer (option) – Top Entry, Side Entry, Bottom Entry with circulation facility.
- Hydraulic lid lifting for fast & easy cleaning of top dish & vessel.
- Load Cell for fast & accurate measuring system.
- Flexible product transfer pipe & fitting in case of a load cell & hydraulic lifting.
- VFD for Anchor agitators, Homogenizer to vary the RPM.
- Specially Design drives assembly.



KEY FEATURES, OUR PRODUCTS AND SERVICES

- Model No. SPEVTD 6 to 96 Tray.
- Stainless steel construction (SS 316L/316/304)
- Easily removable hollow pad type heating shelves
- Dedicated fluid heating and circulating systems
- The thermal efficiency of 60% to 80%
- Provision for Explosion Vent / Rupture Disc on body & Validation Port on Body
- On Chamber Vacuum break valve is provided
- High-quality interlocking for doors to the chamber
- Uniform heating flow
- Branded immersion heaters
- Heavyduty MS / SS door fitted on hinges
- Single and double doors (for large models)
- Complete automation w/ PLC controls (Optional)

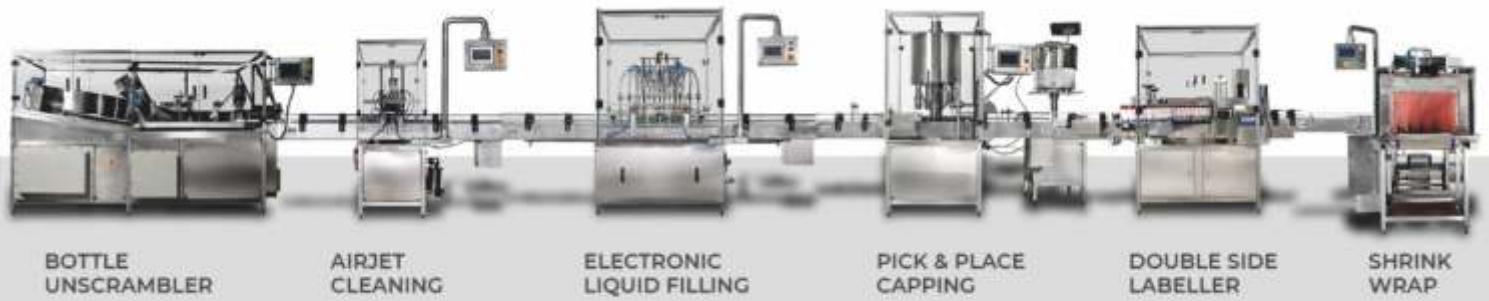

BOTTLE UNSCRAMBLER

- Bottle Unscrambler suitable for any kind of round and flat bottles.
- Output speed from 60 BPM to 300 BPM
- Bottles are discharge upright and properly arranged into existing Conveyor.
- The machine is CRU Design, without any tools, the change of varieties can be finished well.
- PLC Controlled
- Advance container counting & speed detection Sensors with Auto Report Generation.
- Stainless Steel Construction Meets cGMP Design.



AIR JET CLEANING MACHINE

- Bottle Inverter Down & UP position
- No Container – No Air – No Vacuum
- Air Jet Blowing Process In Invert Bottle Position
- Lioniser For Static Removal For Bottles
- Vacuum filter with filter element for Vacuum Blower
- Quick changeover of change parts as per bottles
- Positive air pressure removes the dust from the bottle body.
- Vacuum removes the dust from the bottle
- Stainless Steel Construction Meets cGMP Design

LIQUID FILLING MACHINE

- 4/6/8/12 Filling Station.
- Unique Fill Piston And Servo Nozzle Design Are Driven By Servo System Which Feature Fast And Extremely High Fill Accuracy.
- Various Speed Setting of liquid filling.
- Diving Type and bottom Up Filling Features less Foaming During Filling.
- Simple and Quick Tech Mode Setup.
- Bottle Mouth Holder Design for Special Application.
- PLC Control with Human Machine Interface.
- Error Message Are Instantly Shown on the Touch Panel Which Features Ease Of the Machine Maintenance And Troubleshooting.
- Stainless Steel Construction Meets cGMP Design.



CAPPING MACHINE

- Machine Operating Status Time Setting Monitoring of Abnormalities and Maintenance Instructions Can be Readily Accessed and Easily Operated from the Touch Screen Panel.
- Motor Driven Height Adjustment
- PLC Control With User Friendly HMI
- Capping Mechanism Disengages Automatically When The Present Torque has been reached thus Ensuring proper Capping Without Cap Damage.
- Automatic Cap Elevating System.
- Stainless Steel Construction Meets cGMP Design.



STICKER LABELLING MACHINE

- Self-Set Auto Adjusting Control System
- PLC Control With User Friendly HMI
- Label head with up to $\pm 1\text{mm}$ accuracy (* bottles /containers dependent)
- Product Sensor with Opaque and Clear Bottle Mode
- Variable Speed Separator Wheel
- Up to 200 Bottles per Minute (BPM)
- Stainless Steel Construction Meets cGMP Design



VISUAL INSPECTION MACHINE

- Magnifying Glass with LED Light
- LUX Level > 2000
- Continuous Operation
- Timer based Conveyor
- Variable Frequency Drive
- Stainless Steel Construction Meets cGMP Design



SPARKLER FILTER PRESS



KEY FEATURES, OUR PRODUCTS AND SERVICES

- Design is cGMP with all contact parts AISI 316 L.
- Complete filtration of the batch without any hold up due to reverse flow of liquid.
- The cake, filtered & solids remain enclosed, allowing filtration of toxic, hazardous & explosive liquids.
- No contamination of the internal surface of the shell.
- The top dome structure ensures fast & easy removal of cartridge assembly & cleaning of cake.
- Skid type unit with castor wheel.
- Due to the reverse flow of liquid, there is zero hold up in the filter press that leads to complete filtration of the batch.
- A backwash can provide a facility for continuous operations.
- No scavenger arrangement is required.
- The cake, filter aid & solids remain enclosed, allowing filtration of toxic, hazardous & explosive liquids.
- The top dome structure ensures fast & easy removal of cartridge assembly & cleaning of cake.
- We can increase the number of plates as per toxic levels of input or filtration requirements.
- **AVAILABLE CAPACITIES** : Standard models 8", 14", 18", and 24".

TECHNICAL
 TEAM MATERIALS
 PACKAGING TRUST PROCESS
 DELIVERY
 ISO
 QUALITY
 INSPECTION
 SPEED
 CONTROL
 SUPPLIERS
 QA
 PRODUCTION
 CUSTOMER PARTNERS OPERATIONS
 CHALLENGES ASSEMBLY



SUNRISE PHARMA EQUIPMENTS

📍 **Regd. Office & Factory :** Plot No. 4 Dev Industrial Estate,
 Opp. Pushkar Industrial Estate Phase-1, Vatva GIDC,
 Vatva, Ahmedabad-382 445. Gujarat, India.

☎ +91 98799 97722
 ✉ sales.sunrisepharma@gmail.com
 🌐 www.sunrisepharma.in